

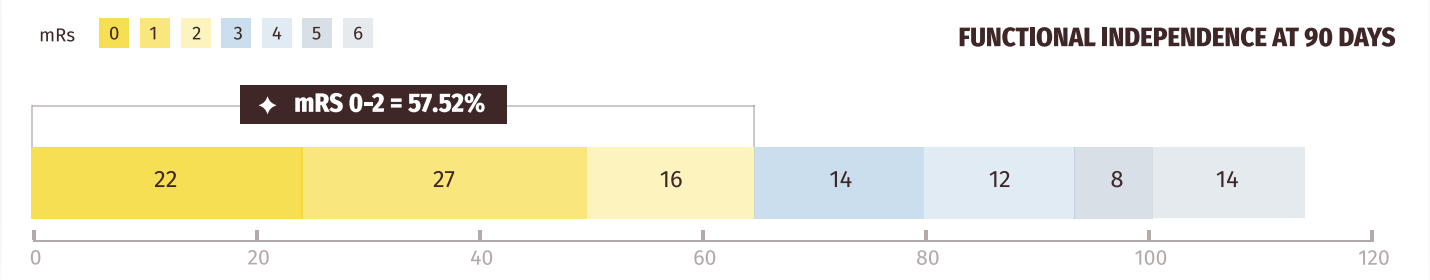
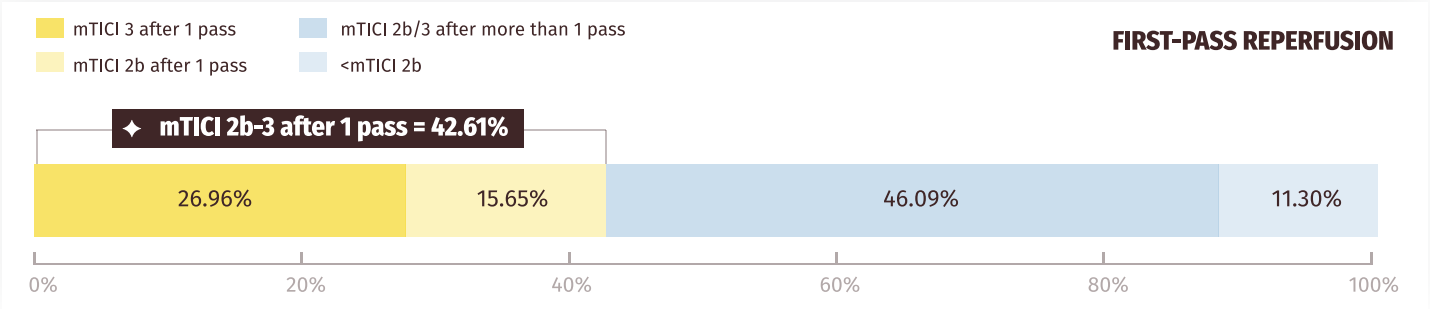


CAPTURE
Key Results

Prospective, Randomized, Controlled, Open Label, Single Blinded,
All China, Multi-Center, Registration Trial of the Thrombectomy
System for the Treatment of Acute Ischemic Stroke
NCT number: NCT04995757.
• 239 patients enrolled and allocated to control arm and NeuroHawk™ arm
• 15 investigate centers enrolled

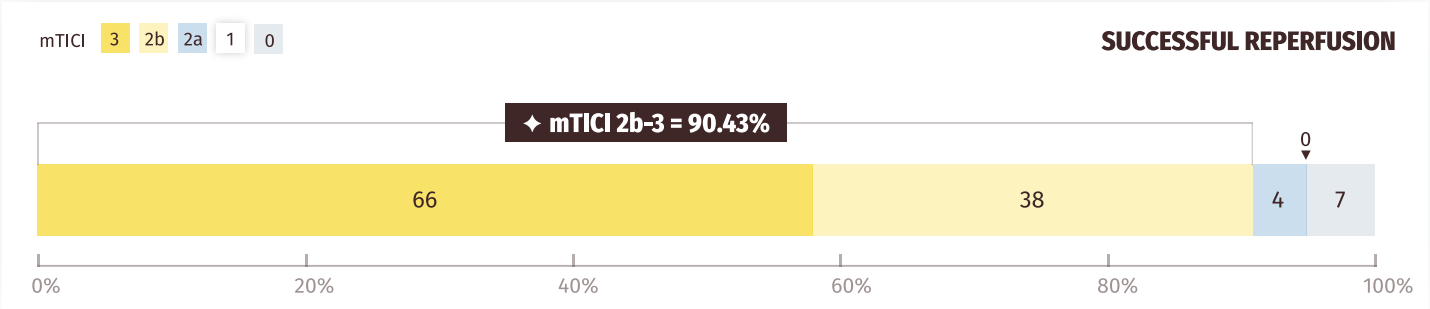


Customized approach, Reliable revascularization

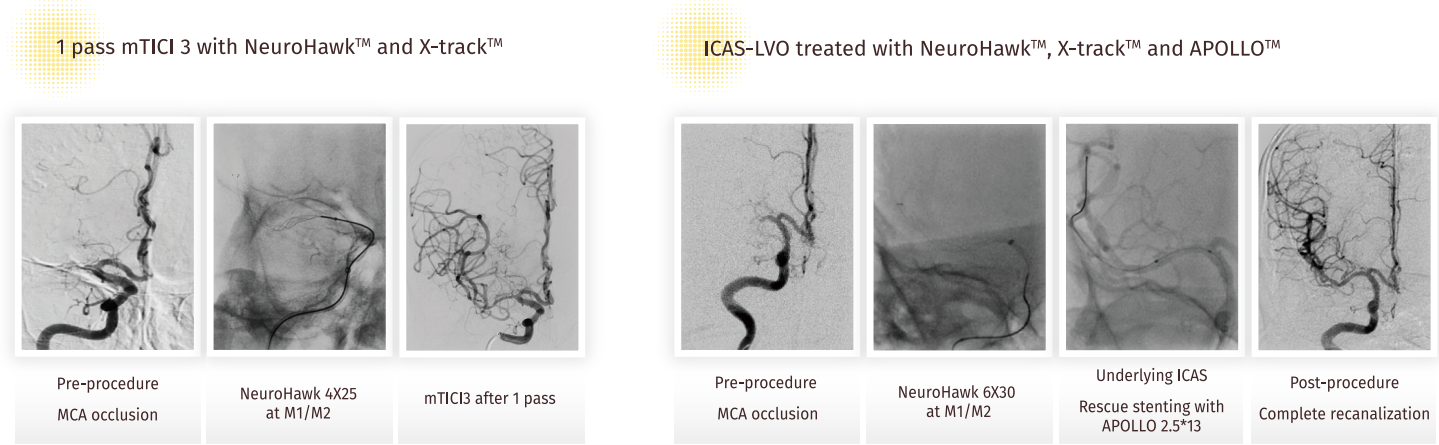

*X-track™ and APOLLO™ are products of MicroPort NeuroTech and may not be available in all regions.

Product Overview

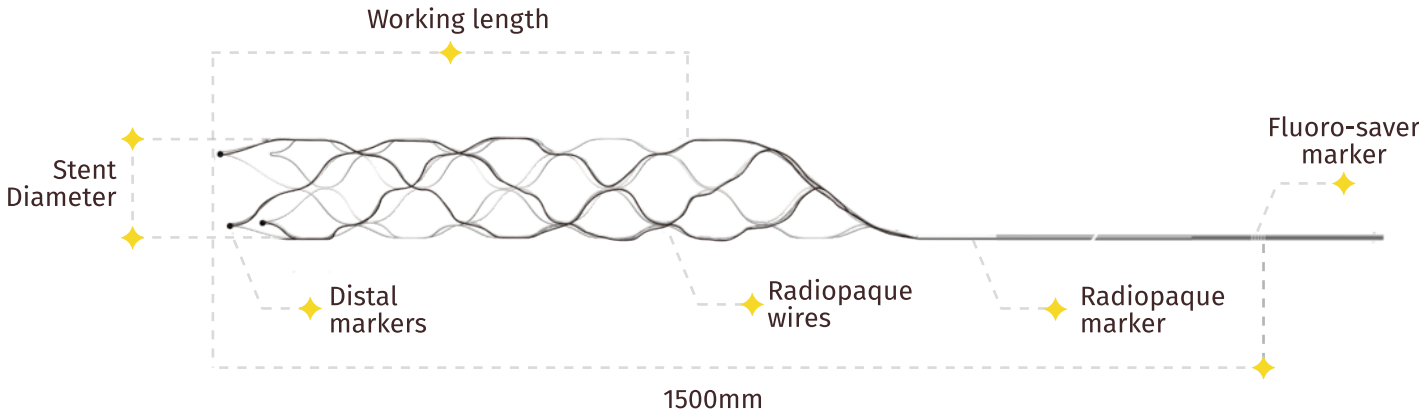


Product Specification

Product REF	Recommended artery diameter	Stent nominal diameter	Working length	Pusher length	Pusher outer diameter	Compatible microcatheter
AIS4025	2-3mm	4.0mm	25mm	200cm	0.018 inch	0.021 inch
AIS6030	3-5mm	6.0mm	30mm			



REDEFINING
THROMBECTOMY
with NEXT generation solutions

MicroPort NeuroTech (Shanghai) Co., Ltd.
1661 Zhangdong Road, Zhangjiang High-tech Park,
Pudong New District, Shanghai, PRC
TEL: (86) (21) 38954600-55200
EMAIL: NT_IBD@microport.com
Website: www.medneurotech.com
LinkedIn: www.linkedin.com/company/microport-neurotech

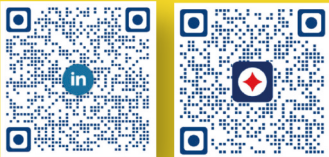


NeuroHawk™
Thrombectomy Device

VERSATILITY

Innovative Hybrid Cell Design

Reliable clot capture and retention of all compositions

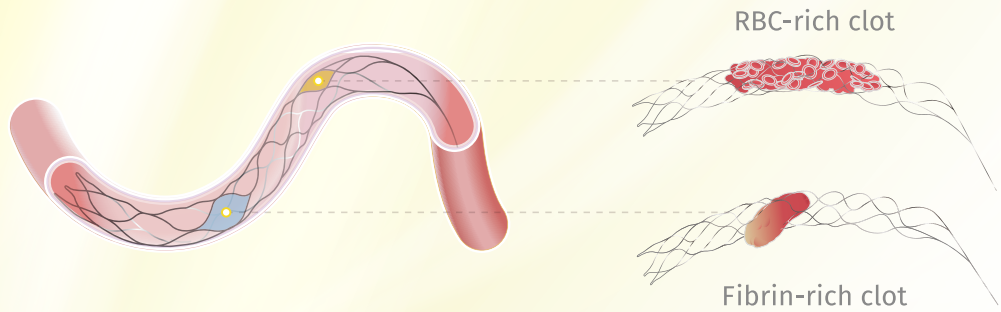


RADIOPACITY

3 Crossing Radiopaque Wires & 3 Distal Markers

Full-length visualization provides


- ◆ Real-time visual feedback of stent expansion and compression
- ◆ Evaluation of clot and lesion types through stent imaging

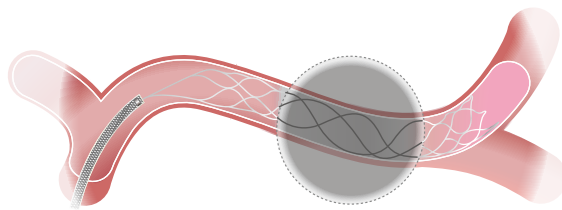


Clot integration
you can **SEE**

DELIVERABILITY

Delivery Force in 0.021" microcatheter
* Data on file.

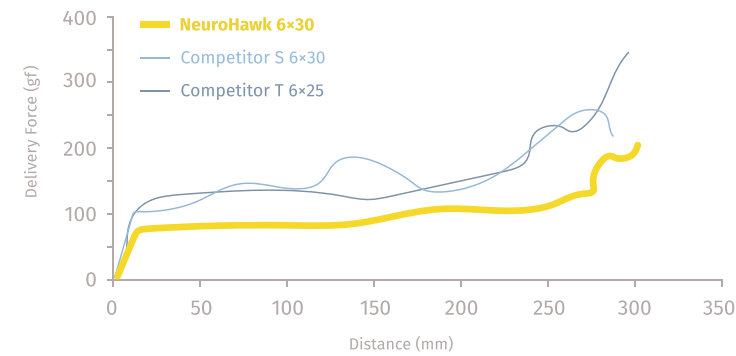

Gradient Core Shaft and Distal Coil Structure
Smooth navigation through tortuous vessels

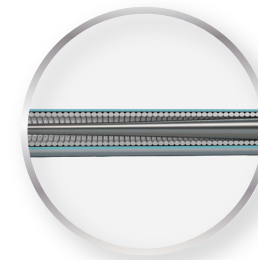


RETRIEVABILITY

Sectional Junction Design

Improved junction flexibility and strength

- ◆ Maximize the effective retrieval force in unfavorable vascular anatomy

Retraction Force(N) from M2 segment

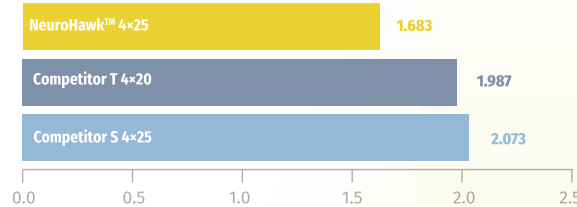

* Data on file.




Clot retrieval
you can **FEEL**

EXPLORE **NEXT** TECHNIQUE

NeuroHawk™

X-track™

NeuroTech offers a comprehensive
stroke portfolio for the treatment
of Acute Ischemic Stroke

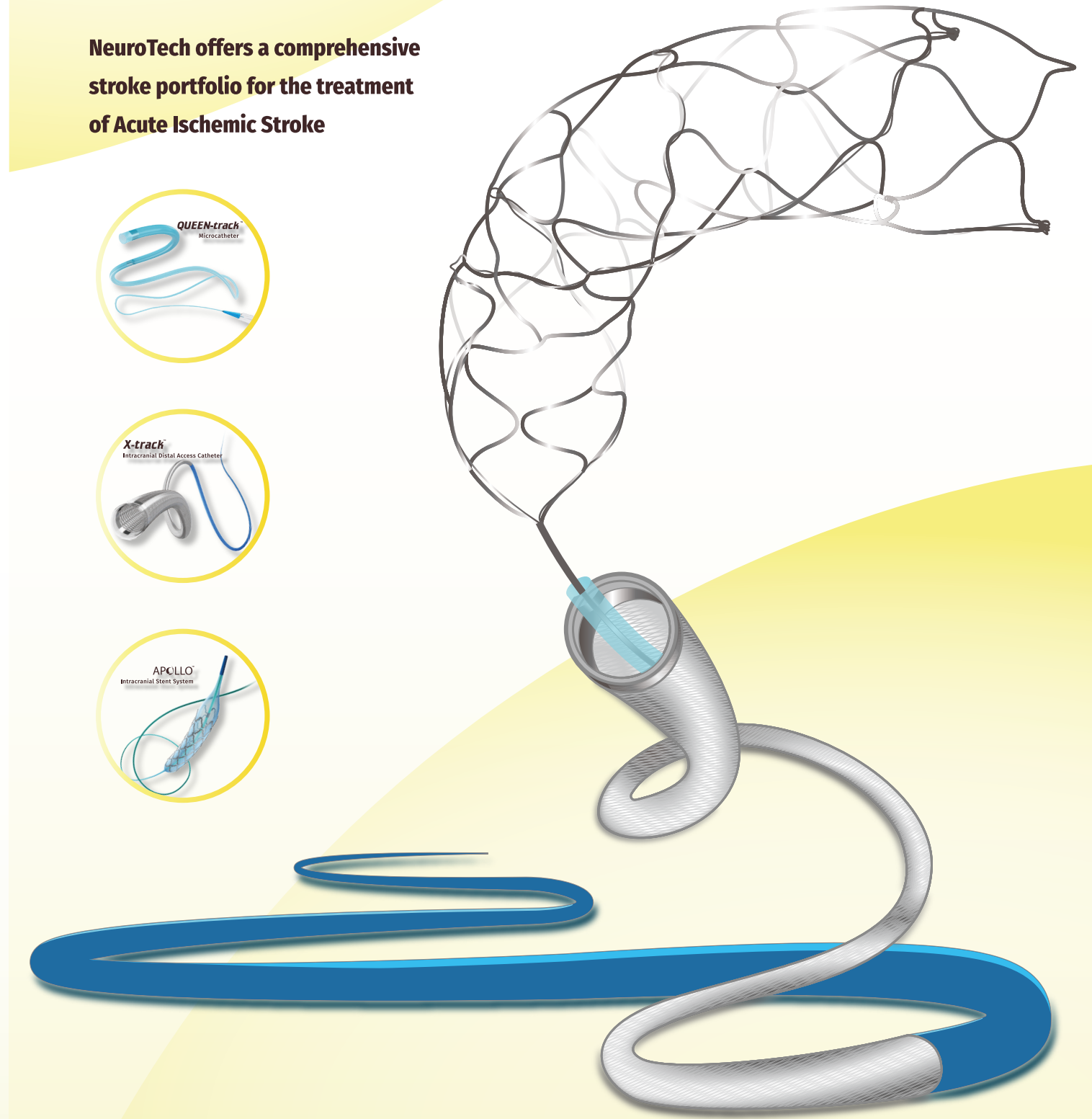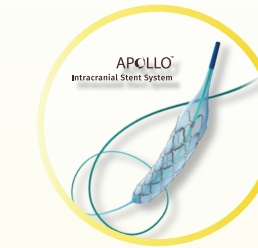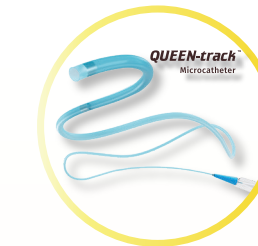

QUEEN-track™, X-track™ and APOLLO™ are products of MicroPort NeuroTech and may not be available in all regions.